


AppliChrom® ABOA SugarSep-Oligo-Na

28/09/2015 Application Note


Oligosaccharidanalytics in Water

AppliChrom ABOA SugarSep-Oligo-Na, 300x8mm Order No SASOLNa103008
"Maltodextrin 19", Eluent: H₂O, 0,25ml/min, 70°C



Total chromatogram

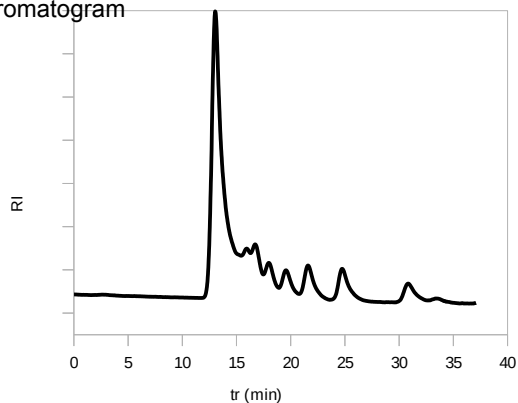



Detail enlargement

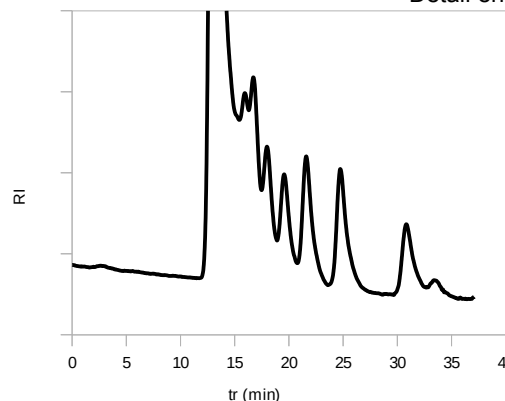

AppliChrom ABOA SugarSep-Oligo-Na, 300x8mm Order No SASOLNa103008
"Maltodextrin 12", Eluent: H₂O, 0,25ml/min, 70°C

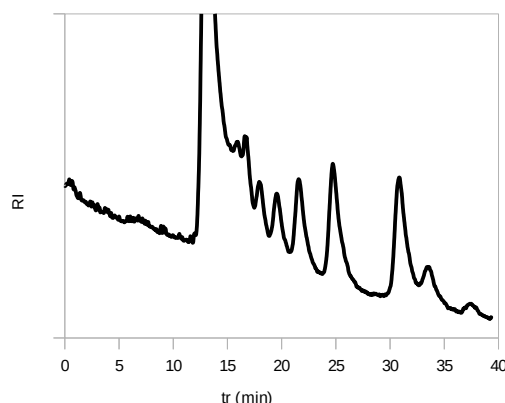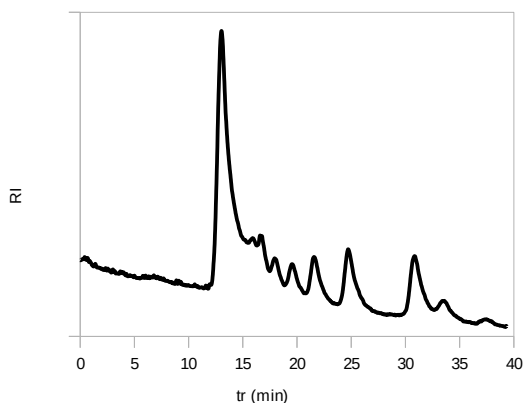

AppliChrom ABOA SugarSep-Oligo-Na, 300x8mm Order No SASOLNa103008
"Maltodextrin 6", Eluent: H₂O, 0,25ml/min, 70°C

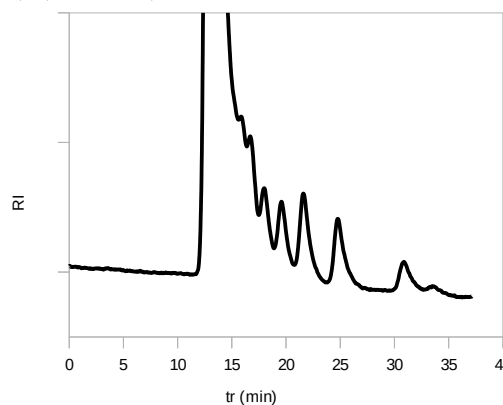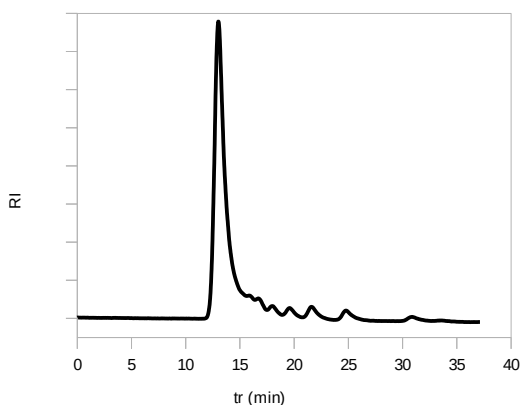

Order Information:

SASOLNa101508	HPLC-Column AppliChrom® ABOA SugarSep-Oligo-Na, 10µ 150mm x 8mm	Price available upon request
SASOLNa103008	HPLC-Column AppliChrom® ABOA SugarSep-Oligo-Na, 10µ 300mm x 8mm	Price available upon request
SASOLNa10508	Precolumn AppliChrom® ABOA SugarSep-Oligo-Na, 10µ, 50mm x 8mm	Price available upon request

AApplication note 09/2015, changes and errors reserved. Not intended for pharmaceutical, diagnostic or medical use, for research and laboratory only. Only terms and conditions of AppliChrom are valid.